
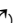


# RASPBERRY PI USER AGENT – FRONTEND MANUAL


## CONNECTING

The front-end can connect to a multitude of Raspberry Pies, either manually by inserting its IP Address into the provided *TextBox* and clicking the *Add* button or by relying on the *RegistryServer* which automatically finds available devices and allows for easily connecting to them. When a *RegistryServer* is available, use the *Refresh*  button for finding available connections.

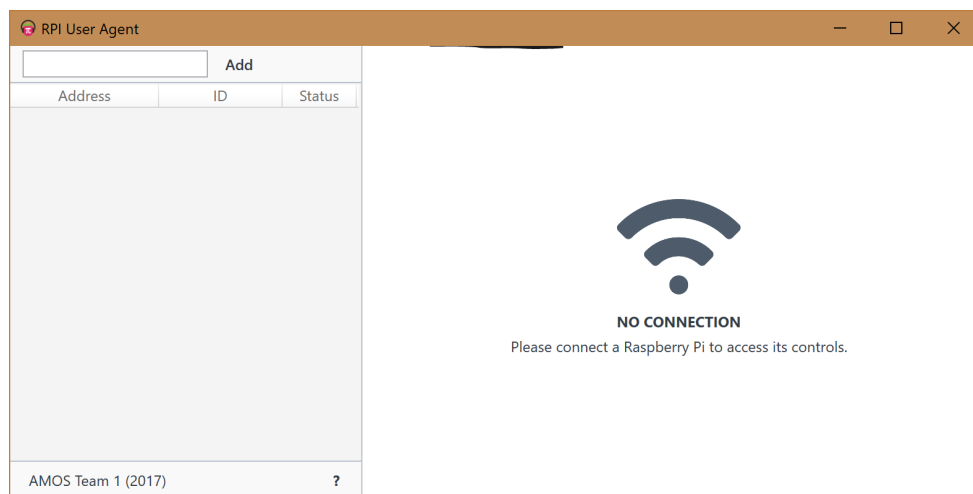

## USER CONTROLS


Once a connection has been successfully established, the *User Controls* for that particular Raspberry will be shown (see below).

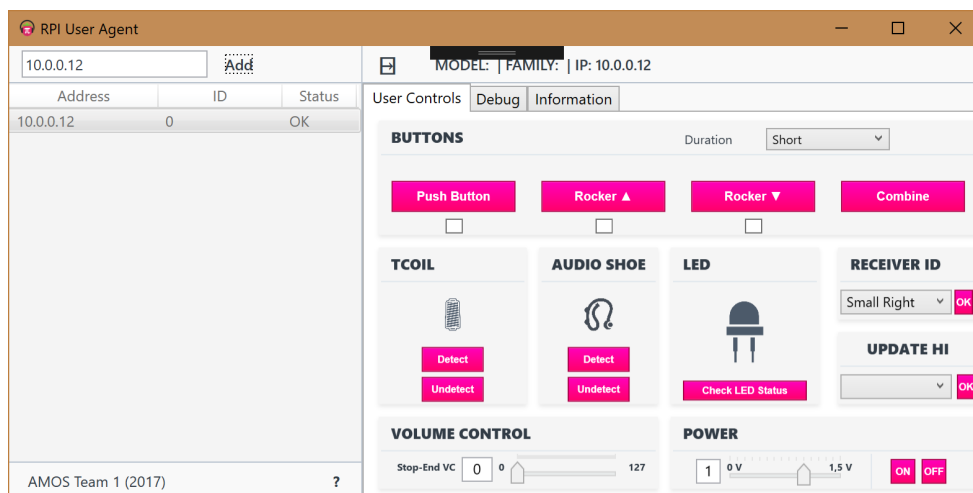

Interacting with this device is easily possible by activating the according hearing aid controls. While most of these are self-explanatory it's worthwhile noting that the buttons can be pushed for different durations. Also, a combination of the *Push Button* with either a *RockerSwitchUp* or *RockerSwitchDown* can be achieved by activating the checkboxes below the buttons and clicking on the *Combine Button*.

Note: It's not possible to press *RockerSwitchUp* and *RockerSwitchDown* simultaneously.

## MINIMAL VIEW

The *User Controls* pane can also be compacted to a minimal appearance where only user controls that are available directly on an actual hearing aid will be shown. This can be activated by pressing the according  button in the menu bar.

